

NON STANDARD EVACUATION TIPS

- **NSNs TO FACILITATE CO-TM LEVEL CASEVAC**

Litter Chains (for M113s): 4010-00-417-2596

Litters(1ea): 6530-01-380-7309



Litter Straps(1ea): 6530-00-784-4205

M113 Red Cross Decals(4ea): 7690-01-103-6312

- **NOTE: Litters and litter straps are class 8 items (medical) and your battalion aid station may have to order these for the company/teams**


M998, UTL VEH, 1.25, 4X4

“The HMMV”

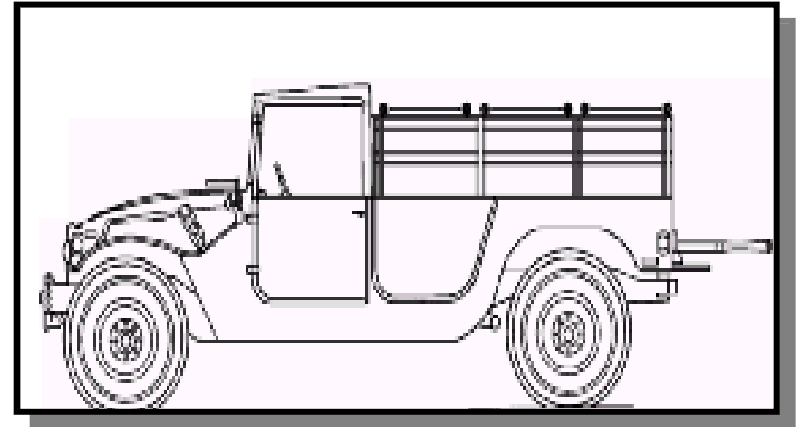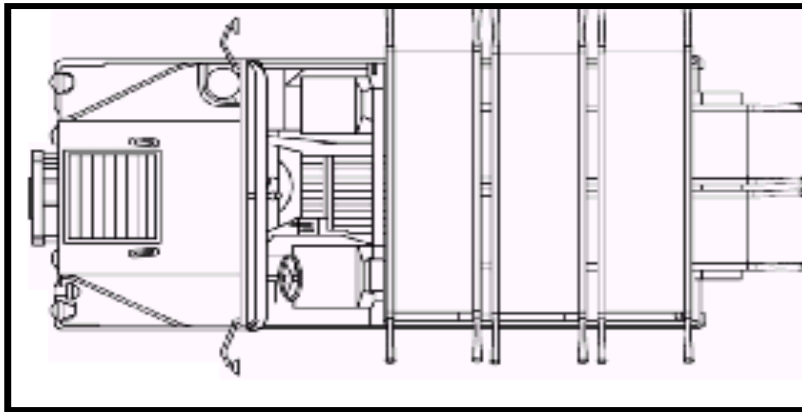


- **PLACE THREE LITTERS SIDE-BY-SIDE ACROSS THE SIDE BOARDS. SECURE THE LITTERS IN PLACE**
- **PLACE TWO LITTERS LENGTHWISE, HEAD FIRST, IN THE BED OF THE TRUCK. SECURE THE LITTERS IN PLACE.**
- **CLOSE THE TAILGATE.**
- **LITTERS ARE UNLOADED IN THE REVERSE ORDER OF LOADING.**

TRUCK, CARGO, 2 1/2 or 5 TON

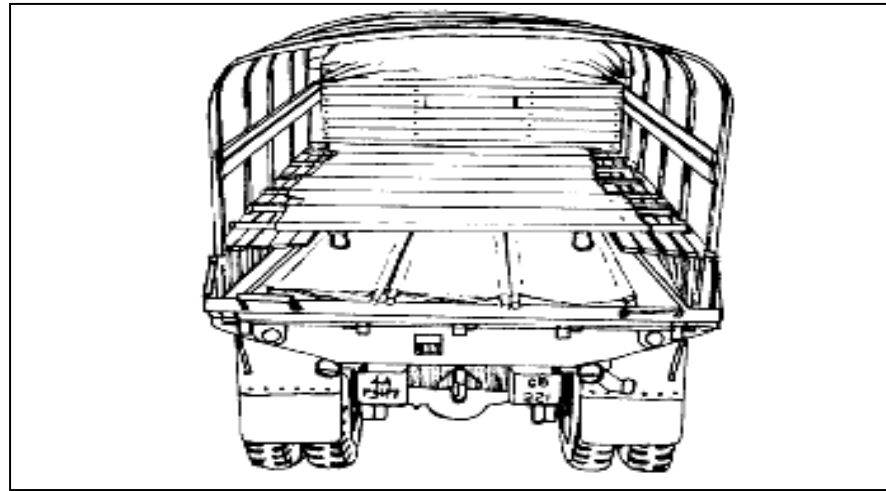


- MAX LITTER CAPACITY IS 12
- LOWER THE SEATS.
- THREE LITTERS CROSSWISE ON THE SEATS AND THREE LITTERS LENGTHWISE ON TRUCK BED (FAR FORWARD AS POSSIBLE)
- SECURE THE LITTERS INDIVIDUALLY ON THE SEATS.
- REPEAT AS REQUIRED
- RAISE THE TAIL GATE


M1095, LMTV

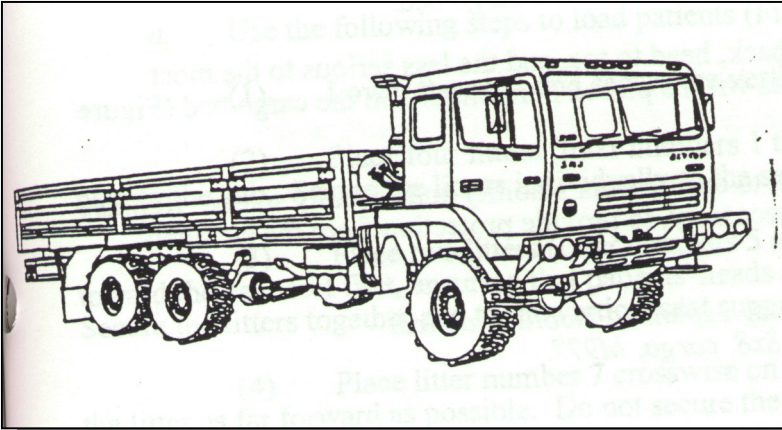


- **MTV LONG WHEEL BASE**
- **12 LITTER CASUALTIES**
- **COMBAT MEDIC RIDES IN THE CENTER**

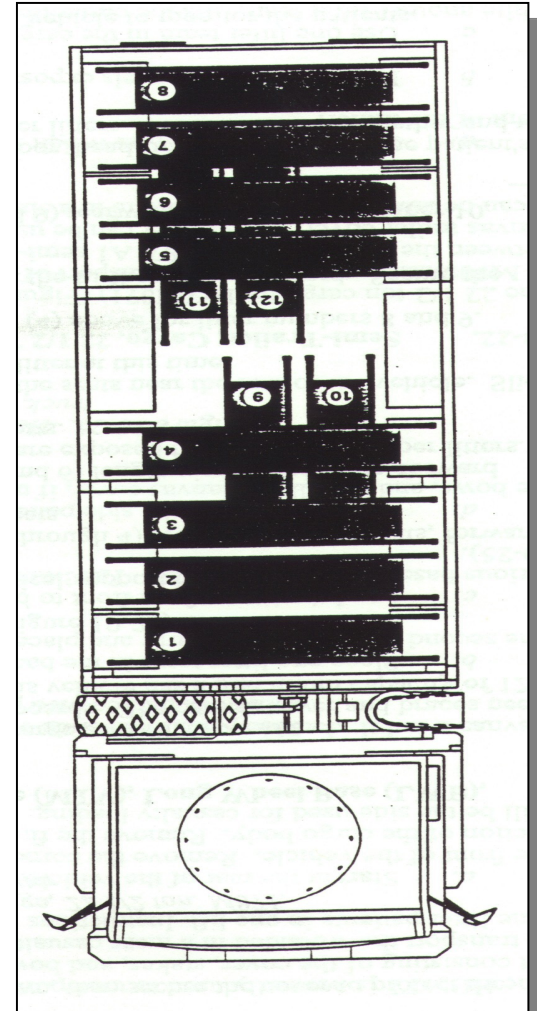



M1093, LMTV

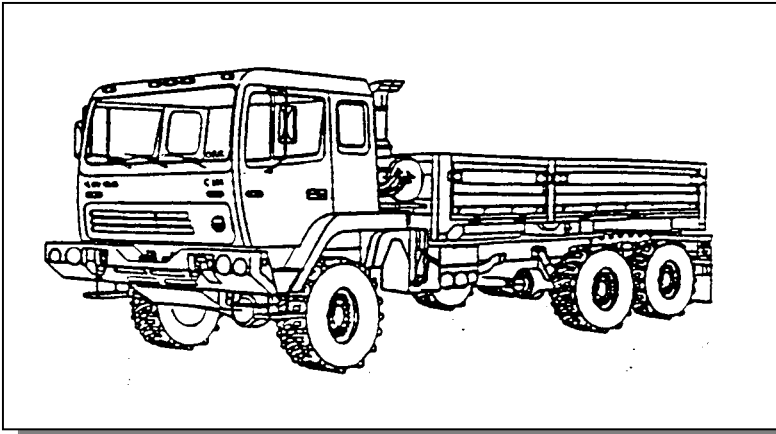

- MTV SHORT WHEEL BASE
- CAN TRANSPORT 8 LITTER
- COMBAT MEDIC IN THE CENTER

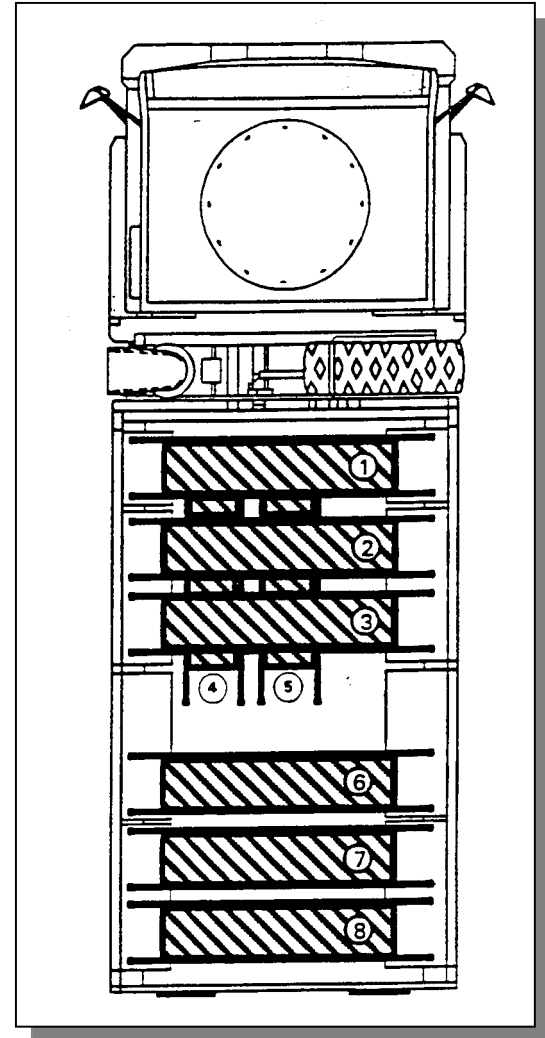

M2, BRADLEY FIGHTING VEHICLE

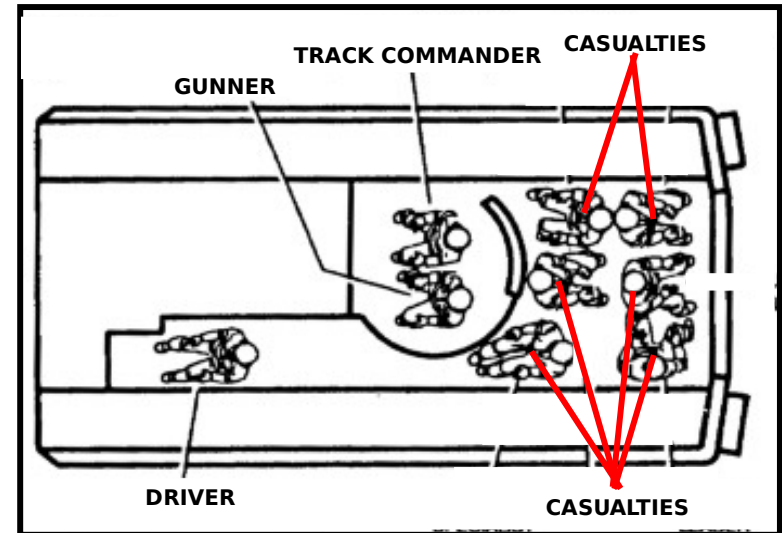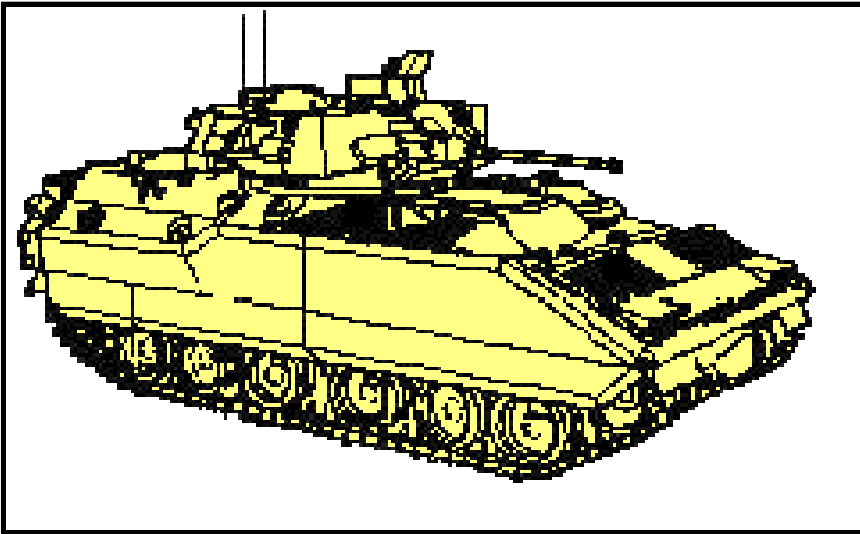

- **BFV CAN TAKE 6 WALKING WOUNDED, NO LITTER CASUALTIES.**
- **THE BFV, QUICKEST AND SAFEST WAY TO EVACUATE CASUALTIES**
- **CASUALTIES SHOULD BE MOVED TO A COVERED AND CONCEALED CCP TO THE REAR OF THE PLATOON**
- **IF ENEMY INDIRECT FIRE PRESENTS THREAT, KEEP CASUALTIES IN THE BFV, TRANSFER TO THE AMBULANCE**